



NaviTip Delivery Tips

When it comes to root canals, Ultradent is with you every step of the way.



NaviTip® Double Sideport

NaviTip® Tip

NaviTip® FX® Tip

Don't Change Your Technique. Make It Easier with NaviTip.

Ultradent's NaviTip tips are the optimal delivery tips for every step in a root canal procedure. Unlike traditional irrigation tips, the cannula of each NaviTip is slightly rigid through the base and center, but flexible at the tip to allow for easy navigation of curved canals.

Simplify Any Procedure

Use NaviTip tips to deliver any endodontic irrigant directly where it's needed—just short of the apex.





NaviTip® 29 ga or 30 ga

- Controlled delivery anywhere in the canal
- Flexible cannula easily navigates curved canals

20pks

50pks



NaviTip® 31 ga Double Sideport Irrigator Tip

- The small cannula navigates the most intricate canals
- Double sideports horizontally express irrigant toward the canal walls, not the apex

20pks

50pks



NaviTip® 30 ga FX® Tips

- Brush tip scrubs, cleans, and irrigates simultaneously
- Flexible cannula easily navigates curved canals
- Listed as an "EXCELLENT" product by a prominent independent research institute*

NaviTip® Tips

	20pk		50pk	
	29 ga	30 ga	29 ga	30 ga
17 mm, White	5112	1249	1378	1421
21 mm, Yellow	5113	1349	1374	1422
25 mm, Blue	5114	1250	1376	1423
27 mm, Green	5115	1354	1377	1424
Assorted	5116	1351	1379	3319

NaviTip® Double Sideport Irrigator Tip

	20pk	50pk
21 mm, Yellow	5121	5122
27 mm, Green	5123	5124

NaviTip® FX® Tips


	20pk
17 mm, White	1452
25 mm, Blue	1454

20pks

NaviTip® adapts to your style

The NaviTip tip delivers low-viscosity materials exactly where they are needed—even in the smallest canals.

The double sideport openings in the NaviTip 31 ga allow your irrigant to be safely expressed against the canal walls rather than toward—or even through—the foramen.

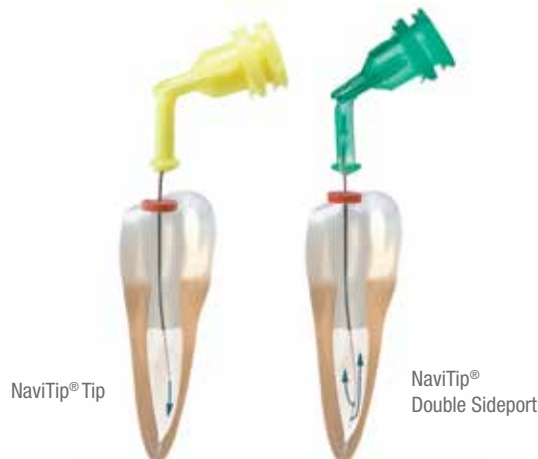

Use NaviTip to deliver products listed below and many other products:

- **ChlorCid® & ChlorCid® V** Sodium Hypochlorite Solutions
- **Ultradent® EDTA** 18% Solution
- **Consepsis® & Consepsis® V** Chlorhexidine Antibacterial Solutions
- **File-Eze®** EDTA Lubricant
- **Ultradent® Citric Acid** 20% Solution
- **UltraCal® XS** Calcium Hydroxide Paste
- **EndoREZ®** Hydrophilic Resin Root Canal Sealer
- **MTAFlow** Repair Cement



*Clinical Research Associates Newsletter. Volume 29, Issue 1, January 2005.